



- With integrated Needle-Free valve
- With PrimeStop minimising exposure and mess
- 20 drops / 1ml
- Roller clamp for flow control
- Spinlock connector for easy kink-free connection
- Suitable for blood and solution administration
- Various configurations for a range of therapy needs.
- Drip chamber filter size: 200µm
- PVC Tubing
- Latex free
- DEHP-free

SANGOFIX®

BLOOD/SOLUTION ADMINISTRATION SETS









SANGOFIX® WITH SAFSITE®

NEEDLE-FREE CAPPED VALVE

Description	Tubing Length	Box Qty	Article No.
 <p>Sangofix® 1 x Safsite®, Y-site with PrimeStop</p>	255 cm	100	4051998
 <p>Sangofix® 1 x Safsite®, 1-way valve with PrimeStop</p>	255 cm	100	4052005
 <p>Sangofix® 1 x Safsite® hand blood pump set with PrimeStop</p>	260 cm	25	4050192
 <p>Sangofix® 2 x Safsite®, 1-way valve, G-clamp, handblood pump set with PrimeStop</p>	260 cm	25	4050193
 <p>Sangofix® 2 x Safsite®, 1-way valve with PrimeStop</p>	255 cm	100	4053710
 <p>Sangofix® 1 x Safsite®, 1-way valve, side arm and G-clamp with PrimeStop</p>	255 cm	100	4052013
 <p>Sangofix® 1 x Safsite® Needle Free Valve</p>	180 cm	100	4052218

SANGOFIX® WITH CARESITE®


NEEDLE-FREE SWABABLE VALVE

Description	Tubing Length	Box Qty	Article No.
 Sangofix® 1 x Caresite, with PrimeStop	255 cm	100	4051998H
 Sangofix® 1 x Caresite, 1-way valve, with PrimeStop	255 cm	100	4052005H
 Sangofix® 2 X Caresite, 2 x 1-way valve, with PrimeStop	255 cm	100	4053710H
 Sangofix® 2 x Caresite, 2 x 1-way valve, hand blood pump set, with PrimeStop	260 cm	25	4050193H
 Sangofix® 1 x Caresite, hand blood pump set, with PrimeStop	260 cm	25	4050192H
 Sangofix® 1-way valve, side arm and G-clamp, with PrimeStop	255 cm	100	4052013H
 Sangofix® Piggyback	150 cm	100	4117301
 Sangofix® 1 x Caresite® Needle Free Valve	180 cm	100	4052218H

Caresite® Features

- Easy-grip barrel design to minimise slip & touch contamination
- Open fluid path design to minimise areas where bacteria can collect
- MRI compatible
- Swabable, smooth surface permits easy cleaning
- Clear design aids flush conformation
- Effortless Luer guide top access
- Caresite® valve integrated into fluid pathway providing continuous flushing

SANGOFIX® PRIMESTOP

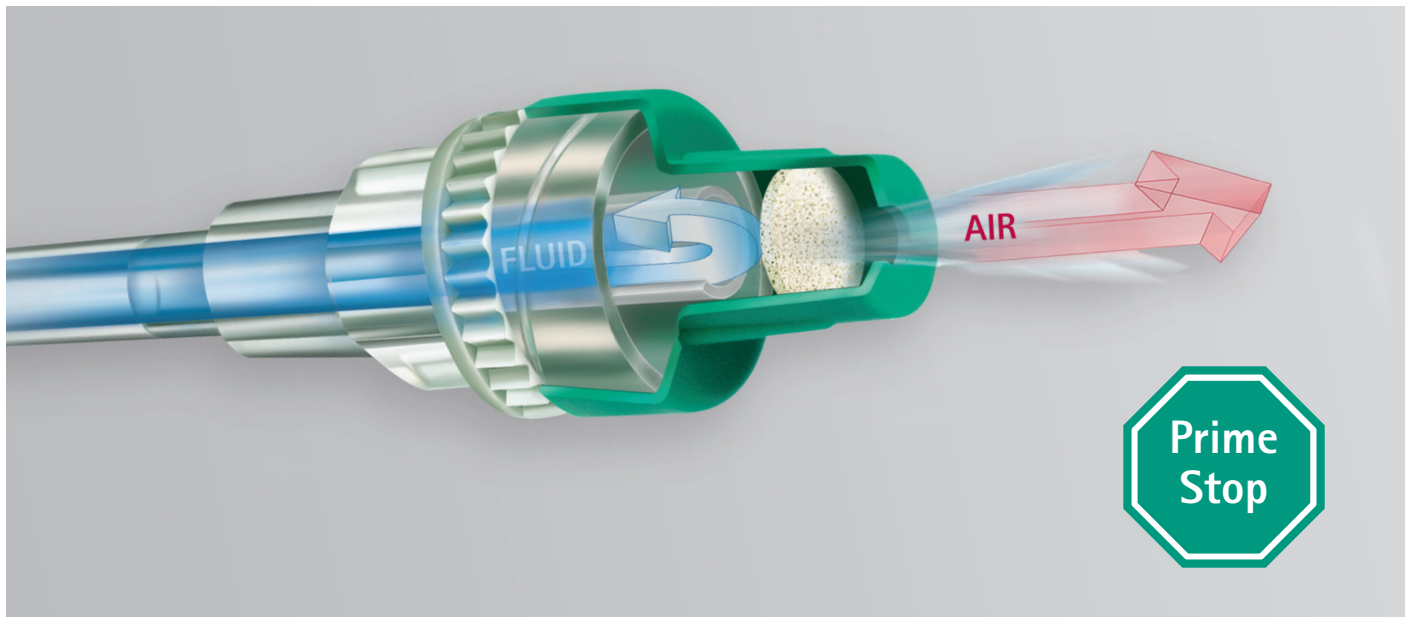

Sangofix® PrimeStop offers cleaner, more efficient priming for IV lines.

A specialised protective cap lined with a hydrophobic membrane stops fluid leaking and protects against contamination during priming.

PrimeStop eliminates the need to expel fluid when priming:

- No more mess
- Reduced contamination and exposure

B. Braun Australia Pty Ltd | Bella Vista, NSW | Tel. 1800 251 705 | info@bbraun.com | www.bbraun.com.au
B. Braun New Zealand | Customer Care 0800 227 286 | Fax (09) 373 5601 | www.bbraun.co.nz