



50 YEARS
YASARGIL CLIPS



NEUROSURGERY

CLIP TRAY

FOR AESCULAP® YASARGIL ANEURYSM CLIPS




CLIP TRAY

FOR AESCULAP® YASARGIL ANEURYSM CLIPS

VALIDATED AND
EASY HANDLING

2 Labor-Medizinprodukte
mit Seriennummer und Barcode

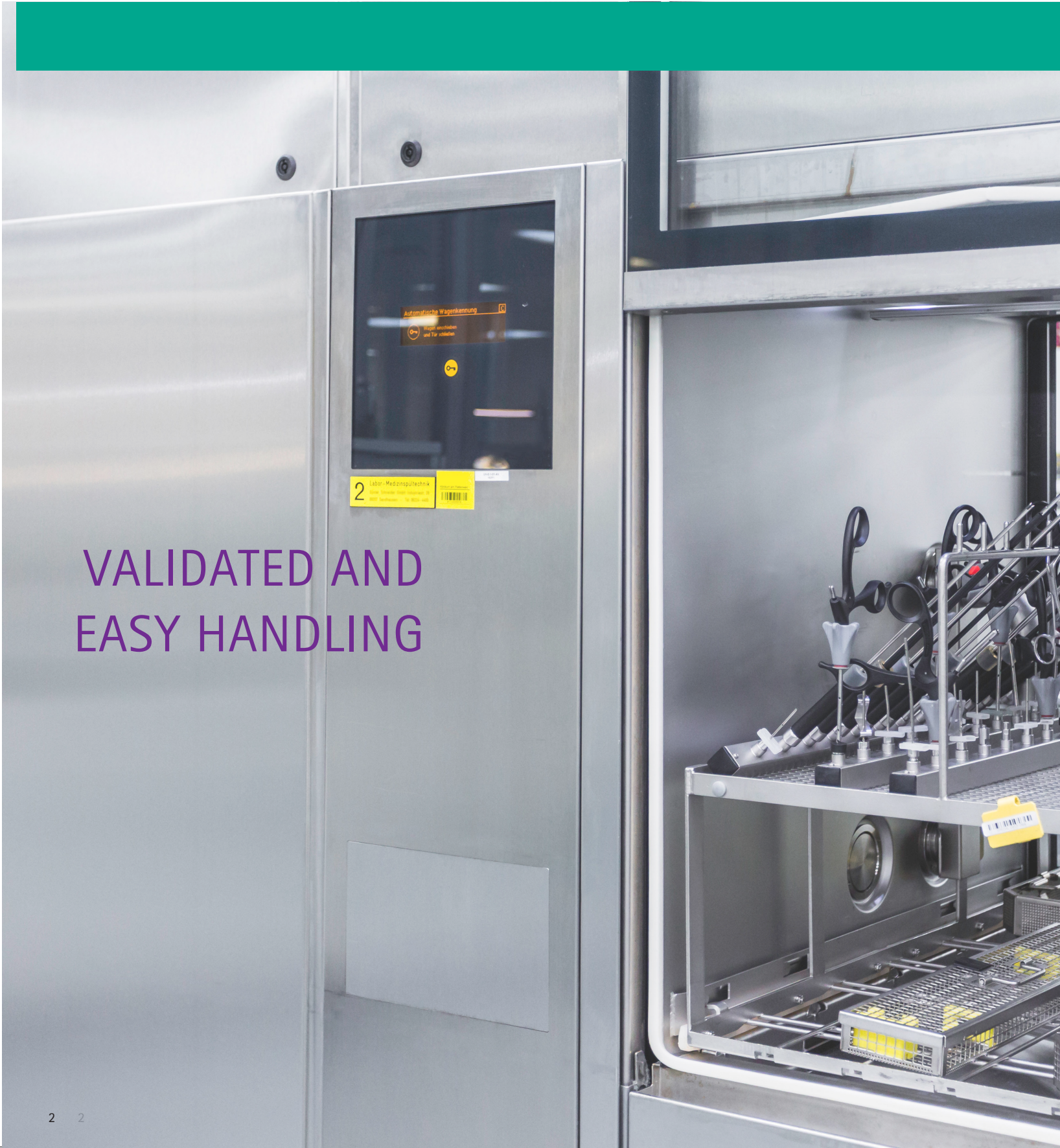



INDIVIDUAL CLIP-ARRANGEMENT

Temporary and permanent clips can be stored in one tray for a rapid retrieval of suitable clips during surgery. The tags can be arranged individually corresponding to your personal preferences.



REDUCED FOOTPRINT TO SAVE SPACE


An increased number of compartments in smaller trays saves space in the OR. The tray is available in two sizes. The small tray offers 21 compartments, the large tray has 42 compartments.

OLD



38 compartments

NEW

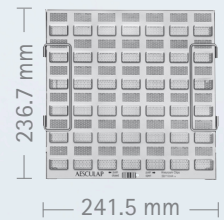

42 compartments

LONG-TERM STABILITY

The tray is made of metal which is more durable than current plastic trays even after hundreds of reprocessing cycles.

TRACEABILITY OF EACH CLIP

A Trace Plate makes the tracking of each clip easier. This Trace Plate is used to store the patient labels of the permanent clips that are in the clip tray.

CLIP TRAY

FOR AESCULAP® YASARGIL ANEURYSM CLIPS




VALIDATED & EASY HANDLING OF THE CLIPS

The tray provides a validated, clear and easy way to organize, store and **reprocess** your AESCULAP® YASARGIL Aneurysm clips.



VALIDATED REPROCESSABILITY OF ANEURYSM CLIPS INSIDE THE TRAY

- 1. Validated reprocessing process:** Cleaning and sterilization of temporary and permanent aneurysm clips inside the tray is validated. Tests show good results with regard to cleaning, biocompatibility and endotoxines.
- 2. Good drying properties:** The tray is made of stainless steel that prevents the appearance of residual water.
- 3. Improved perforation design:** The special design of the perforations of the tray offers a good water permeability.
- 4. Long-term stability:** Usage of high-alkaline cleanser is possible.

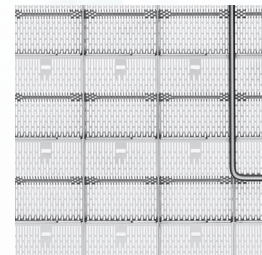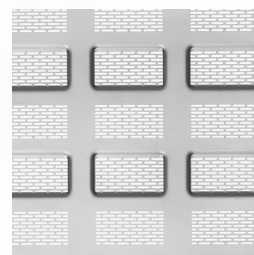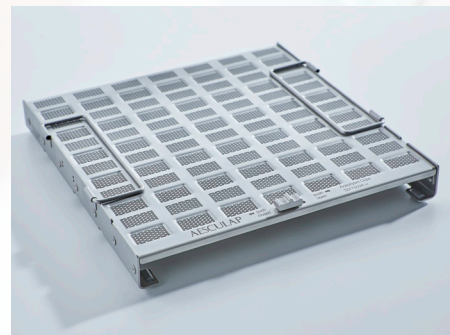

IDENTIFICATION TAGS



For your individual clip arrangement, a corresponding clip tag is available for each clip. The tag can be attached to the tray easily. The identification tags contain the article number, the blade length, the spring size as well as a picture in original size of each clip.


For an easy allocation of the clips to the tags, the identification tags have the same color coding as the AESCULAP® YASARGIL Aneurysm clips. Tags for temporary clips are yellow, tags for permanent clips are grey.

LONGER COMPARTMENTS FOR LONG CLIPS



The AESCULAP® "long clips" require longer compartments. This is being achieved by removing a separator and replacing it by a shorter separator.

TRACE PLATE FOR PATIENT LABELS

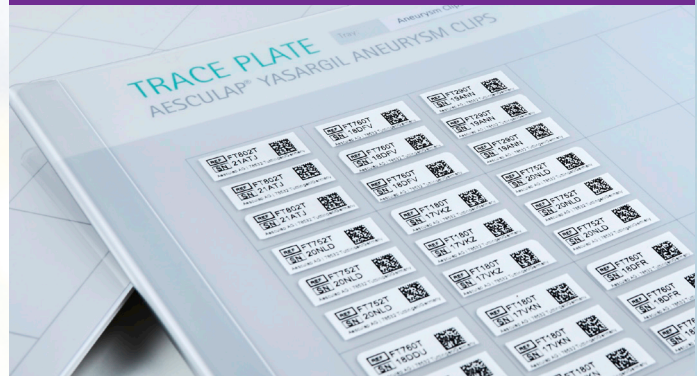

A Trace Plate for the storage of the patient labels is available. It represents the design of the small / large tray.

IDENTIFICATION LABELS





The identification labels make an individual marking of the trays possible.


CLIP TRAY

PRODUCT INFORMATION

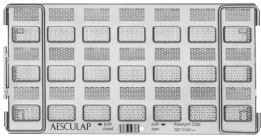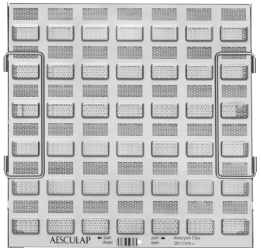

CLIP TRAY

Art. No.	Image	Description	Size
FT008R		Aneurysm Clip Tray small	241.5 x 125 x 37.8 mm
FT009R		Aneurysm Clip Tray large	241.5 x 236.7 x 37.8 mm



IDENTIFICATION TAGS

Art. No.	Image	Description
Identification Tags Kennzeichnungsschilder		The article numbers of the identification tags are made up of the clip article number plus the suffix "-ID". Example: FT750T-ID


TRACE PLATES

Art. No.	Image	Description
C02902		Trace Plate for Aneurysm Clip Tray small (FT008R)
C03002		Trace Plate for Aneurysm Clip Tray large (FT009R)

SEPARATORS

Art. No.	Image	Description	Note
FT008205		Separator short for Aneurysm Clip Tray small (FT008R)	Separators are used to create longer compartments for long clips.
FT009205		Separator for Aneurysm Clip Tray large (FT009R)	

IDENTIFICATION LABELS FOR TRAY

Art. No.	Image	Description
Identification labels for tray		Identification labels are available in different colours and with / without lettering. For more information see catalog C40402.



B. Braun Australia Pty Ltd | Bella Vista, NSW | Tel.1800 251 705 | info@bbraun.com | www.bbraun.com.au
B. Braun New Zealand | Customer Care 0800 227 286 | Fax (09) 373 5601 | www.bbraun.co.nz

BAUS AE E0595 05/20 | ©2020 B. Braun Australia